**NÁVOD K POUŽITÍ
NÁVOD NA POUŽITIE
INSTRUKCJA OBSŁUGI
INSTRUCTION MANUAL
HASZNÁLATI ÚTMUTATÓ**

**PEP 817
Trézlidé**

Gallet

**Osobní váha s tělesnou analýzou
Osobná váha s telesnou analýzou
Waga łazienkowa z analizą ciała
Bathroom body analysis scale
Személyi mérleg test analízissal**

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

HOUSEHOLD USE ONLY – Tylko do zastosowania w gospodarstwie domowym.
DO NOT IMMERSE IN WATER OR OTHER LIQUIDS – Nie zanurzać do wody lub innych cieczy. TO AVOID DANGER OF SUFFOCATION, KEEP THIS PLASTIC BAG AWAY FROM BABIES AND CHILDREN. DO NOT USE THIS BAG IN CRIBS, BEDS, CARRIAGES OR PLAYPENS. THIS BAG IS NOT A TOY.

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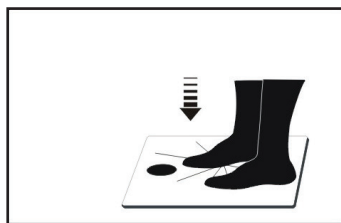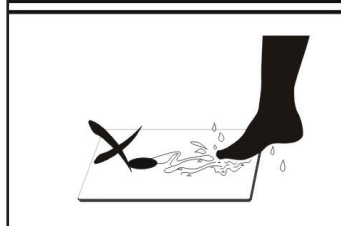Dear customer, thank you for purchasing our product. Please read the operating instructions carefully before putting the appliance into operation and keep these instructions including the warranty, the receipt and, if possible, the box with the internal packing.

GENERAL SAFETY INSTRUCTIONS


- The appliance is designed exclusively for private use and for the envisaged purpose.
- This appliance is not fit for commercial use.
- Keep it away from sources of heat, direct sunlight, humidity (never dip it into any liquid) and sharp edges. Do not use the appliance with wet hands.
- The appliance can use by people (including children) whose physical, sensual or mental inability or lack of experience and knowledge prevents safe use of the appliance if they are supervised.
- Never use the appliance if it malfunctions, if it fell to the ground and got damaged or if it fell into water. In these cases bring the appliance to a specialised service to have the safety and functions of the device checked.
- Keep in mind that the data read from this scale are approximate and inaccurate in comparison to results of actual medicinal analyses. The exact determination of biometrical characteristics can be done only by a doctor via medical methods.
- Do not use during pregnancy, with symptom of edema and during dialysis treatment.
- Do not use the fat analysis feature on this scale if you are using one of the following devices: pace maker, electrocardiograph and/or other medical device implanted in the body or device used as support.
- If you should leave the workplace, always switch the machine off.
- Do not try to repair the appliance on your own. Always contact an authorized technician.
- Pay careful attention to the following "Special Safety Instructions for this unit".

SPECIAL SAFETY INSTRUCTIONS FOR THIS UNIT

- Operate the appliance only in accordance with the specifications on the rating plate.
- Do not put the appliance on the wet places.
- Clean the scale with a slightly damp cloth. Do not immerse the scale in the water or use chemical / abrasive agents.
- Always use the scale on a hard, flat surface. DO NOT use on carpet
- DO NOT use the scale on wet and slippery floor. DO NOT step on the edge of the scale. The scale platform is slippery when wet. Keep it dry!
- Always weigh yourself with the same scales places on the same surface. Uneven floors may affect the reading.
- Weigh yourself without clothing and footwear, before meals and always at the same time of day.
- To make sure accuracy of testing, please DO NOT use the scale in environment with strong electro-magnetic field.
- DO NOT use the scale on wet and slippery floor. DO NOT step on the edge of the scale. The scale platform is slippery when wet. Keep it dry!
- DO NOT use this scale on a place with intense shake. Keep stand still on scale when testing.
- Do not drop the device, or jump on the device when using. Avoid strong impact and sudden shake as it is a precision measuring device. The construction is delicate and sensitive to all unwanted external impact.
- All text in other languages and pictures shown on the packaging, or the product itself, are translated and explained at the end of this language mutation.
- The producer is not responsible for damage caused by improper use of the appliance (e.g. breaking the board, etc.) and its guarantee does not cover the appliance in the case of failure to comply with the safety warnings above.

DESCRIPTION OF THE CONTROLS

- P1 LCD Display
 P2 Electrode
 P3 SET key
 P4 Down key ▼
 P5 UP key ▲

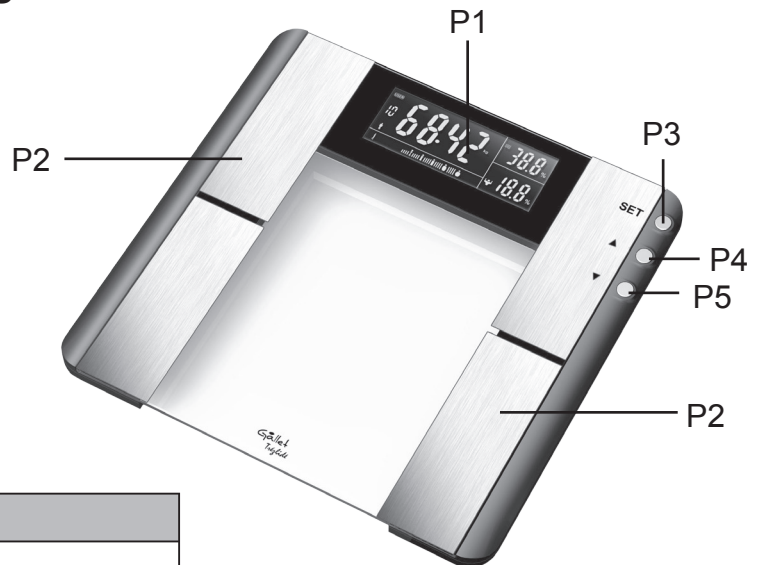Display symbol	Meaning
kg (lb/st)	Body weight
USER	User memory space
	Male
	Female
cm (Ft:in)	Body height
Age	Age
	Unathletic type
	Average type
	Athletic Type
	Body fat display
	Body water display
	Muscle mass display
	Bone weight display
%	Percentage – body fat / body water / muscle mass / bone weight
	Analysis bar: Extremely underweight – under - weight – normal weight – over - weight – obese
KCAL	Daily kilocalorie requirement
Increase	Weight gain in kg (lb/st) compared with previous measurement
Dncrease	Weight loss in kg (lb/st) compared with previous measurement


General Instructions

Gallet Glass Body scale applies BIA (Bio-impedance Analysis) technology. A small amount of weak current flows through the human body so as to detect the bio-impedance and estimate body fat, body water, muscle mass and bone mass. The electrical current is small and may not be felt. This BIA technology is cheap, safe, non-invasive, toxic-free and harmless. It also possesses the characteristics of simple operation and abundant information.


The current mentioned above is less than 1mA. However, please be aware that anyone with an wearable or implantable medical electronic instrument, such as a pacemaker, must avoid using this device. The intended use of this device is for adult's indoor use only.

Insert the Batteries

- Open the battery door in the back of the scale.
- Insert the batteries (2 x 1.5V AAA) into the battery compartment according to the polarity indications marked inside the compartment.
 - * The digits "88888" will be shown on the LCD.

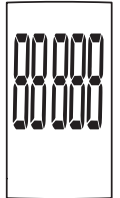- Close the battery door and wait until the digits "0.0kg" are shown on the LCD.

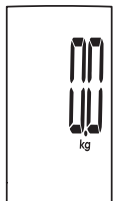CAUTION

- When the symbol **LO** appears, the device will power off in four seconds. Then you shall replace with a new set of batteries. Please replace all four batteries in the same time. Do NOT mix the old batteries with the new one.
- Worn batteries are hazardous waste. Do NOT dispose of them together with the household garbage. Please refer to the local ordinances and recycling instructions regarding disposal of the worn batteries.

Set Up Your Profile

The glass body scale supports multiple users (Up to 10). You may follow below instructions to assign User ID and set up your own profile, including Gender, Stature, and Age.

1. Assigning User ID

- With batteries correctly installed, press "**SET**" key to enter setting.
- The system will request User ID selection first. "**USER**" blinks.
 - The operator may press the function key ▼ or ▲ to select User ID among 0 - 9.
- Press "**SET**" key to confirm User ID.

2. Setting Gender

- After confirming User ID, the system will divert to Gender setting automatically.
- The portrait  blinks. The operator may press the function key ▼ or ▲ to select Gender and Status.
- Press "**SET**" key to confirm User ID.

3. Setting Age


- After confirming Gender, the system will divert to Age setting automatically.
- The digits "**AGE**" blinks. The operator may press the function key ▼ or ▲ to increase or decrease the numeral. (Range: 10 to 85 years old)
- You may press and hold the function key ▼ or ▲ for fast changing the numeral.
- Press "**SET**" key to confirm Age.

4. Setting Stature

- After confirming Age, the system will divert to Stature setting automatically.
- The digits of units of height blinks. The operator may press the function key ▼ or ▲ to increase or decrease the numeral. (Range: 75 cm to 225 cm).
- You may press and hold the function key ▼ or ▲ for fast changing the numeral.
- Press “**SET**” key to confirm stature.

5. Setting type

The default symbol for the “unathletic” type flashes. Which type do you consider yourself to be? Select your type using the ▼ or ▲ button:

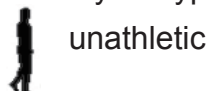unathletic

aver-age

athletic

Confirm your selection with SET.

When the LCD displays “0.0kg”, you may start measuring.

Select Measurement Unit

To change the settings to “lb” or “st” as well as “Ft:in”, proceed as follows:

- Switch the scales on by applying gentle pressure on the glass surface with your foot or toe.
- Press the ▲ button, until the required unit of weight is displayed. The unit for the height also automatically changes.

First Measurement

- Please process the user setting before your first measurement.

(Refer to Set Up Your Profile for more details)

STEP 1: Step on the platform barefooted or press **SET**.


STEP 2: The preset stored position data appears on the display (USER flashes). Use the ▼ or ▲ button to select the stored position at which your data is saved.

STEP 3: Wait for the USER display and 0. 0 kg (0.0 lb/st) to stop flashing.

STEP 4: Stand on the contact surfaces of the scales with both feet.

STEP 5: After a few seconds, the following results will be shown:


- your weight / the body fat content in % / the body water content in % / the results analysis in the form of a bar
- the calorie requirement per day / the muscle mass content in % / the bone mass content in %
- The weight gain (INCREASE) / the weight loss (DECREASE) in kg (lb/st) in comparison to the previous measurement / The preset gender / The preset type / The preset height /- The preset age (AGE)

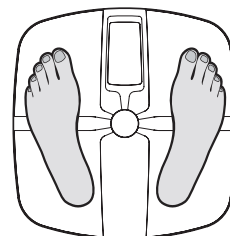Afterwards the determined data will be shown repeated. Then the scale will automatically turn off.

- If it fails to identify the possible User ID, the LCD will only display the weight data.
(To find out the solutions, please refer to Troubleshooting for more details).

Body Fat Indicator

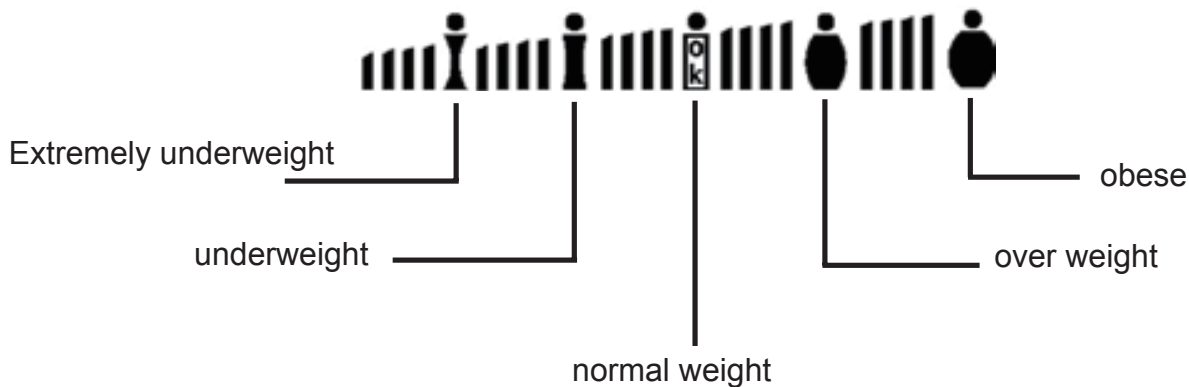Display error Messages

Display symbol	Meaning
8888	Self test (when switching on) to check that the numeric display functions properly
— — — —	The scales stand faultily
Lo	Battery is low
Err	The percentage of body fat content lies above or below the limits of the measuring scale
EEEE	The weight lies above the allowed maximum weight.

When Measuring ...

Problem	Root Cause	Solution
Abnormal measuring results: - Too high; OR - Too low; OR - Huge difference between two recent measurement.	Incorrect posture	Please step on the platform barefooted and stand still.
	The device is located on the soft ground such as a carpet OR on a rugged surface.	Please place the device on a flat, hard surface.
	Cold body that may results in bad blood circulation.	Warm up your hands and feet to resume blood circulation and then measure again.
	Cold Electrodes.	Place the device in a warm room for a while and then measure again.
	Either your hands or your feet are too dry.	Wipe your feet with a damp cloth, keeping them slightly damp when starting measurement.

Problem	Root Cause	Solution
No display on LCD when the device powers on.	Batteries not yet installed.	Install the batteries. (Please refer to Insert the Batteries)
	Worn batteries.	Replace all four batteries in the same time. Please purchase the authorized batteries for replacement.
CANNOT proceed to analyze body fat, total body water, muscle mass and bone mass.	Step onto the platform wearing socks or shoes.	Please keep barefooted during the measurement, and keep full contact with the electrodes as well.
	The system cannot identify the possible User ID with most similar data.	Please assign a User ID following the instruction in Set Up Your Profile.
	The user fails to select the User ID from what the system found.	Please assign a User ID following the instruction in Set Up Your Profile.
The device powers off automatically.	Low battery.	Replace all four batteries in the same time. Please purchase the authorized batteries for replacement.

Table of Body Fat Level (Unit: %)

Age	Women		Man		
	Body Fat	Body Water	Body Fat	Body Water	
<= 30	4,0 - 20,5	66,0 - 5,7	4,0 - 15,5	66,0 - 28,1	Underweight
	20,6 - 25	54,6 - 51,6	15,6 - 20,0	58,0 - 55,0	Normal
	25,1 - 30,5	51,5 - 47,8	20,1 - 24,5	54,9 - 51,9	Overweight
	30,6 - 50,0	47,7 - 34,4	24,6 - 50,0	51,8 - 34,4	Obese
> = 30	4,0 - 25,0	66,0 - 51,6	4,0 - 19,5	66,0 - 55,4	Underweight
	25,1 - 30,0	51,5 - 48,1	19,6 - 24,0	55,3 - 52,3	Normal
	30,1 - 35,0	48,0 - 44,7	24,1 - 28,5	52,2 - 49,2	Overweight
	35,1 - 50,0	44,6 - 34,4	28,6 - 50,0	49,1 - 34,4	Obese

CLEANING AND MAINTENANCE


- 1) Use a dry, soft fabric to clean the device. Or a lightly damp soft fabric with mild detergent to wipe clean the device. Do not use strong chemical solution such as alcohol, benzene for cleaning.
- 2) It could be slippery when the unit platform surface is wet, or semi-damp. Ensure the weighing platform surface is dry before using for safety reason.

TECHNICAL SPECIFICATIONS

- Bathroom body analysis scale
- Analysis of body fat, body water, + muscle mass, bone mass, recommended KCAL
- Measurement Range: 2,5 kg - 150 kg / st: 0,394 - 23:6 / 5,512 - 330,7 lb
- Measurement accuracy 100 g
- Measurement units 00 kg/lb/st
- Large LCD display
- LCD display with white backlit of symbols
- Memory for 10 person
- Weighs immediately - just step on the scales - it will automatically turn on, ready to weigh
- Auto off after 30 sec
- Overload indicator, low battery indicator
- Operation on 2 pieces of AAA batteries

We reserve the right to change technical specifications.

Old electrical appliances, used batteries and accumulators disposal

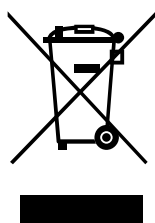This symbol appearing on the product, on the product accessories or on the product packing means that the product must not be disposed as household waste. When the product/ battery durability is over, please, deliver the product or battery (if it is enclosed) to the respective collection point, where the electrical appliances or batteries will be recycled. The places, where the used electrical appliances are collected, exist in the European Union and in other European countries as well. By proper disposal of the product you can prevent possible negative impact on environment and human health, which might otherwise occur as a consequence of improper manipulation with the product or battery/ accumulator. Recycling of materials contributes to protection of natural resources. Therefore, please, do not throw the old electrical appliances and batteries/ accumulators in the household waste. Information, where it is possible to leave the old electrical appliances for free, is provided at your local authority, at the store where you have bought the product. Information, where you can leave the batteries and accumulators for free, is provided to you at the store, at your local authority.

AUTO STEP ON – just step on the scale and stand still, the scale will automatically turn ON and start taking measurements.

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